


BioVendor

Research
and Diagnostic Products

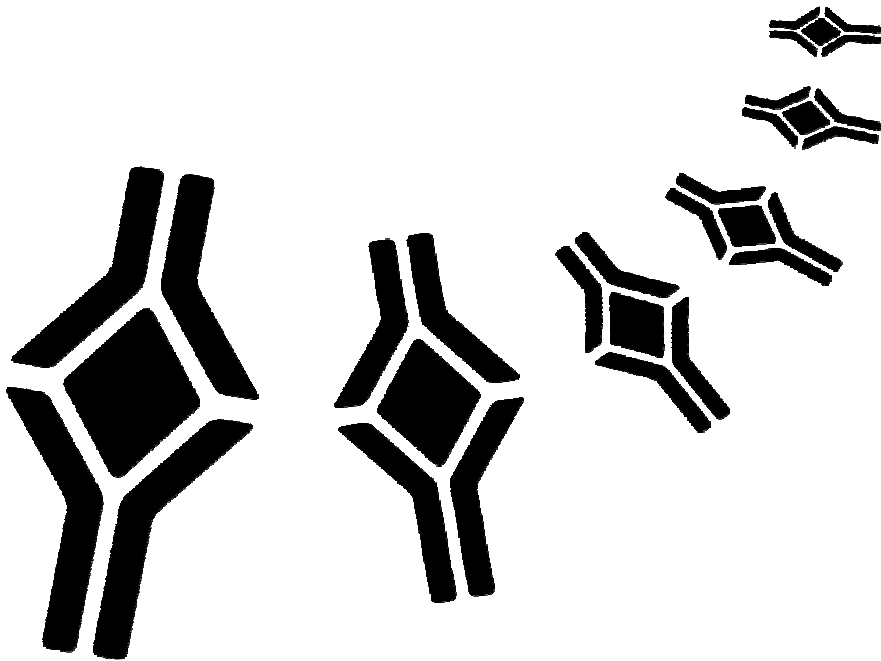

17 α -HYDROXYPROGESTERONE (17 α -OHP) ELISA

Product Data Sheet

Cat. No.: RCAN-P-320R

For Research Use Only

CONTENTS

1.	INTENDED USE	3
2.	PRINCIPLE OF THE TEST	3
3.	CLINICAL APPLICATIONS	3
4.	PROCEDURAL CAUTIONS AND WARNINGS	3
5.	LIMITATIONS	4
6.	SAFETY CAUTIONS AND WARNINGS	4
6.1	Potential Biohazardous Material	4
7.	SPECIMEN COLLECTION AND STORAGE	5
8.	SPECIMEN PRETREATMENT	5
9.	REAGENTS AND EQUIPMENT NEEDED BUT NOT PROVIDED	5
10.	REAGENTS PROVIDED	5
11	CALCULATIONS	8
12	TYPICAL TABULATED DATA	8
13	TYPICAL CALIBRATOR CURVE	8
14	PERFORMANCE CHARACTERISTICS	9
14.1	SENSITIVITY	9
14.2	SPECIFICITY (CROSS REACTIVITY)	9
15.	EXPECTED NORMAL VALUES	11
16.	REFERENCES	11

➤➤ This kit is manufactured by:
BioVendor – Laboratorní medicína, a.s.

➤➤ Use only the current version of Product Data Sheet enclosed with the kit!

1. INTENDED USE

For the direct quantitative determination of 17 α -hydroxyprogesterone (17 α -OHP) in human serum by an enzyme immunoassay. For *in vitro* use only. Multispecies specificity.

2. PRINCIPLE OF THE TEST

The principle of the following enzyme immunoassay test follows the typical competitive binding scenario. Competition occurs between an unlabeled antigen (present in standards, control and patient samples) and an enzyme-labelled antigen (conjugate) for a limited number of antibody binding sites on the microwell plate. The washing and decanting procedures remove unbound materials. After the washing step, the enzyme substrate is added. The enzymatic reaction is terminated by addition of the stopping solution. The absorbance is measured on a microtiter plate reader. The intensity of the colour formed is inversely proportional to the concentration of 17 α -OHP in the sample. A set of standards is used to plot a standard curve from which the amount of 17 α -OHP in patient samples and controls can be directly read.

3. CLINICAL APPLICATIONS

The steroid 17 α -hydroxyprogesterone is produced by the adrenal cortex and gonads. 17 α -OHP has little progestational activity, but it is of intense clinical interest because it is the immediate precursor to 11-desoxycortisol (Cpd-S). Cpd-S is produced by the 21-hydroxylation of 17 α -OHP. Measurement of 17 α -OHP is a useful indirect indicator of 21-hydroxylase activity. In congenital 21-hydroxylase deficiency, the most common variety of congenital adrenal hyperplasia (CAH), 17 α -OHP is secreted in abundant excess. Measurement of 17 α -OHP is therefore valuable in the initial diagnosis of CAH.

4. PROCEDURAL CAUTIONS AND WARNINGS

1. Users should have a thorough understanding of this protocol for the successful use of this kit. Reliable performance will only be attained by strict and careful adherence to the instructions provided.
2. Control materials or serum pools should be included in every run at a high and low level for assessing the reliability of results.
3. When the use of water is specified for dilution or reconstitution, use deionized or distilled water.
4. In order to reduce exposure to potentially harmful substances, gloves should be worn when handling kit reagents and human specimens.
5. All kit reagents and specimens should be brought to room temperature and mixed gently but thoroughly before use. Avoid repeated freezing and thawing of reagents and specimens.
6. A calibrator curve must be established for every run.
7. The control should be included in every run and fall within established confidence limits.
8. Improper procedural techniques, imprecise pipetting, incomplete washing as well as improper reagent storage may be indicated when assay values for the control do not reflect established ranges.
9. When reading the microplate, the presence of bubbles in the microwells will affect the optical densities (ODs). Carefully remove any bubbles before performing the reading step.

10. The substrate solution (TMB) is sensitive to light and should remain colourless if properly stored. Instability or contamination may be indicated by the development of a blue colour, in which case it should not be used.
11. When dispensing the substrate and stopping solution, do not use pipettes in which these liquids will come into contact with any metal parts.
12. To prevent contamination of reagents, use a new disposable pipette tip for dispensing each reagent, sample, standard and control.
13. Do not mix various lot numbers of kit components within a test and do not use any component beyond the expiration date printed on the label.
14. Kit reagents must be regarded as hazardous waste and disposed of according to national regulations.

5. LIMITATIONS

1. All the reagents within the kit are calibrated for the direct determination of 17 α -OHP in human serum. The kit is not calibrated for the determination of 17 α -OHP in saliva, plasma or other specimens of human or animal origin.
2. Do not use grossly hemolyzed, grossly lipemic, icteric or improperly stored serum.
3. Any samples or control sera containing azide or thimerosal are not compatible with this kit, as they may lead to false results.
4. Only calibrator A may be used to dilute any high serum samples. The use of any other reagent may lead to false results.
5. The results obtained with this kit should never be used as the sole basis for a clinical diagnosis. For example, the occurrence of heterophilic antibodies in patients regularly exposed to animals or animal products has the potential of causing interferences in immunological tests. Consequently, the clinical diagnosis should include all aspects of a patient's background including the frequency of exposure to animals/products if false results are suspected.

6. SAFETY CAUTIONS AND WARNINGS

6.1 Potential Biohazardous Material

Human serum that may be used in the preparation of the standards and control has been tested and found to be non-reactive for Hepatitis B surface antigen and has also been tested for the presence of antibodies to HCV and Human Immunodeficiency Virus (HIV) and found to be negative. However no test method can offer complete assurance that HIV, HCV and Hepatitis B virus or any infectious agents are absent. The reagents should be considered a potential biohazard and handled with the same precautions as applied to any blood specimen.

6.2 Chemical Hazards

Avoid contact with reagents containing TMB, hydrogen peroxide and sulfuric acid. If contacted with any of these reagents, wash with plenty of water. TMB is a suspected carcinogen.

7. SPECIMEN COLLECTION AND STORAGE

Approximately 0.1 mL of serum is required per duplicate determination. Collect 4-5 mL of blood into an appropriately labelled tube and allow it to clot. Centrifuge and carefully remove the serum layer. Store at 4°C for up to 24 hours or at -10°C or lower if the analyses are to be done at a later date. Consider all human specimens as possible biohazardous materials and take appropriate precautions when handling.

8. SPECIMEN PRETREATMENT

This assay is a direct system; no specimen pretreatment is necessary.

9. REAGENTS AND EQUIPMENT NEEDED BUT NOT PROVIDED

1. Precision pipettes to dispense 25, 50, 100, 150 and 300 µL
2. Disposable pipette tips
3. Distilled or deionized water
4. Plate shaker
5. Microwell plate reader with a filter set at 450nm and an upper OD limit of 3.0 or greater* (see assay procedure step 10).

10. REAGENTS PROVIDED

10.1 Rabbit Anti-17 α -OHP Antibody Coated Microwell Plate-Break Apart Wells- Ready To Use.

Contents: One 96 well (12x8) polyclonal antibody-coated microwell plate in a resealable pouch with desiccant.

Storage: Refrigerate at 2-8°C

Stability: 12 months or as indicated on label.

10.2 17 α -OHP -Horseradish Peroxidase (HRP) Conjugate Concentrate – X50

Contents: 17 α -OHP-HRP conjugate in a protein-based buffer with a non-mercury preservative.

Volume: 0.3 mL/vial

Storage: Refrigerate at 2-8°C

Stability: 12 months or as indicated on label.

Preparation: Dilute 1:100 in assay buffer before use (eg. 20 µL of HRP in 2 mL of assay buffer). If the whole plate is to be used dilute 180 µL of HRP in 18 mL of assay buffer. Discard any that is left over.

10.3 17 α -OHP Calibrators - Ready To Use.

Contents: Seven vials containing 17 α -OHP in a protein-based buffer with a non-mercury preservative. Prepared by spiking buffer with a defined quantity of 17 α -OHP.

*Listed below are approximate concentrations, please refer to vial labels for exact concentrations.

Calibrator	Concentration	Volume/Vial
Calibrator A	0 ng/mL	2.0 mL
Calibrator B	0.15 ng/mL	0.6 mL
Calibrator C	0.5 ng/mL	0.6 mL
Calibrator D	1.5 ng/mL	0.6 mL
Calibrator E	3 ng/mL	0.6 mL
Calibrator F	7.5 ng/mL	0.6 mL
Calibrator G	20 ng/mL	0.6 mL

Storage: Refrigerate at 2-8°C

Stability: 12 months in unopened vials or as indicated on label. Once opened, the standards should be used within 14 days or aliquoted and stored frozen. Avoid multiple freezing and thawing cycles.

10.4 Control - Ready To Use.

Contents: Two vials containing 17 α -OHP in a protein-based buffer with a non-mercury preservative. Prepared by spiking buffer with a defined quantity of 17 α -OHP. Refer to vial label for expected value and acceptable range.

Volume: 0.6 mL/vial

Storage: Refrigerate at 2-8°C

Stability: 12 months in unopened vial or as indicated on label. Once opened, the control should be used within 14 days or aliquoted and stored frozen. Avoid multiple freezing and thawing cycles.

10.5 Wash Buffer Concentrate - X10

Contents: One bottle containing buffer with a non-ionic detergent and a non-mercury preservative.

Volume: 50 mL/bottle

Storage: Refrigerate at 2-8°C

Stability: 12 months or as indicated on label.

Preparation: Dilute 1:10 in distilled or deionized water before use. If the whole plate is to be used dilute 50 mL of the wash buffer concentrate in 450 mL of water.

10.6 Assay Buffer - Ready To Use*.

Contents: One bottle containing a protein-based buffer with a non-mercury preservative.

Volume: 20 mL/bottle

Storage: Refrigerate at 2-8°C

Stability: 12 months or as indicated on label.

10.7 TMB Substrate - Ready To Use.

Contents: One bottle containing tetramethylbenzidine and hydrogen peroxide in a non-DMF or DMSO containing buffer.

Volume: 16 mL/bottle

Storage: Refrigerate at 2-8°C

Stability: 12 months or as indicated on label.

10.8 Stopping Solution - Ready To Use.

Contents: One bottle containing 1M sulfuric acid.

Volume: 6 mL/bottle

Storage: Refrigerate at 2-8°C

Stability: 12 months or as indicated on label.

ASSAY PROCEDURE

Specimen Pretreatment: **None.**

All reagents must reach room temperature before use. Calibrators, controls and specimen samples should be assayed in duplicate. Once the procedure has been started, all steps should be completed without interruption.

1. Prepare working solution of the 17 α -OHP-HRP conjugate and wash buffer.
2. Remove the required number of microwell strips. Reseal the bag and return any unused strips to the refrigerator.
3. Pipette 50 μ L of each calibrator, control and specimen sample into correspondingly labelled wells in duplicate.
4. Pipette 150 μ L of the 17 α -OHP-HRP conjugate working solution into each well (We recommend using a multichannel pipette).
5. Incubate on a plate shaker (approximately 200 rpm) for 1 hour at room temperature.
6. Wash the wells 3 times with 300 μ L of diluted wash buffer per well and tap the plate firmly against absorbent paper to ensure that it is dry (The use of a washer is recommended).
7. Pipette 150 μ L of TMB substrate into each well at timed intervals.
8. Incubate on a plate shaker for 10-15 minutes at room temperature (or until calibrator A attains dark blue colour for desired OD).
9. Pipette 50 μ L of stopping solution into each well at the same timed intervals as in step 7.
10. Read the plate on a microwell plate reader at 450nm within 20 minutes after addition of the stopping solution.

11 CALCULATIONS


- 1 Calculate the mean optical density of each calibrator duplicate.
- 2 Draw a calibrator curve on semi-log paper with the mean optical densities on the Y-axis and the calibrator concentrations on the X-axis. If immunoassay software is being used, a 4-parameter curve is recommended.
- 3 Calculate the mean optical density of each unknown duplicate.
- 4 Read the values of the unknowns directly off the calibrator curve.
- 5 If a sample reads more than 20 ng/mL then dilute it with calibrator A at a dilution of no more than 1:8. The result obtained should be multiplied by the dilution factor.

12 TYPICAL TABULATED DATA

Calibrator	OD 1	OD 2	Mean OD	Value (ng/mL)
A	2.531	2.558	2.544	0
B	2.357	2.326	2.341	0.15
C	1.776	1.749	1.762	0.5
D	0.957	1.016	0.986	1.5
E	0.668	0.732	0.700	3
F	0.348	0.360	0.354	7.5
G	0.198	0.201	0.199	20
Unknown	0.823	0.825	0.824	2.15

13 TYPICAL CALIBRATOR CURVE

Sample curve only. **Do not** use to calculate results.

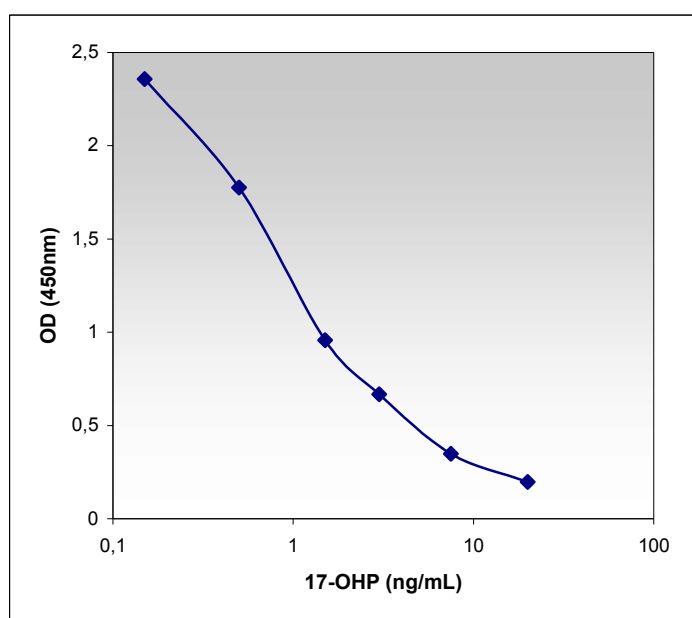

14 PERFORMANCE CHARACTERISTICS

14.1 SENSITIVITY

The lower detection limit is calculated from the standard curve by determining the resulting concentration of the mean OD of calibrator A (based on 10 replicate analyses) minus 2 SD. Therefore, the sensitivity of the 17 α -OHP ELISA kit is

0.11 ng/mL.

14.2 SPECIFICITY (CROSS REACTIVITY)

The following compounds were tested for cross-reactivity with the 17 α -OHP ELISA kit with 17 α -OHP cross-reacting at 100%.

Steroid	%Cross-Reactivity
17 α -Hydroxyprogesterone	100
Progesterone	0.25
11-Desoxycortisol	0.77
DHEA	<0.1
Corticosterone	<0.1
Cholesterol	<0.1
Pregnenolone	<0.1
Pregnenolone-SO ₄	<0.1
Prednisone	<0.1

14.3 INTRA-ASSAY PRECISION

Three samples were assayed 20 times, on the same calibrator curve. The results (in ng/mL) are tabulated below:

Sample	Mean	SD	CV%
1	1.91	0.096	5.0
2	5.68	0.287	5.1
3	10.48	0.602	5.7

14.4 INTER-ASSAY PRECISION

Three samples were assayed 10 times over a period of four weeks. The results (in ng/mL) are tabulated below:

Sample	Mean	SD	CV%
1	1.69	0.116	6.8
2	5.69	0.487	8.6
3	9.38	0.464	4.9

14.5 RECOVERY

Spiked samples were prepared by adding defined amounts of 17 α -OHP to three patient serum samples. The results (in ng/mL) are tabulated below:

Sample	Obs.Result	Exp.Result	Recovery%
1	1.82	-	-
0.15*	2.61	2.75	95
0.4*	4.66	4.8	97
1.6*	13.81	14.8	93
2	2.79	-	-
0.15*	3.57	3.57	100
0.4*	5.90	5.66	104
1.6*	15.47	15.66	99
3	3.70	-	-
0.15*	4.29	4.3	100
0.4*	5.942	6.4	93
1.6*	16.49	16.4	101

* - Amount added (ng/0.1 mL)

14.6 LINEARITY

Three patient serum samples were diluted with calibrator A. The results (in ng/mL) are tabulated below:

Sample	Obs.Result	Recovery%
1	3.52	-
1:2	1.50	86
1:4	0.88	101
1:8	0.46	107
2	5.47	-
1:2	2.7	99
1:4	1.21	89
1:8	0.66	97
3	18.29	-
1:2	8.75	96
1:4	4.54	99
1:8	2.27	100


15. EXPECTED NORMAL VALUES

As for all clinical assays each laboratory should collect data and establish their own range of expected normal values.

Group	Range (ng/mL)
Newborn	
5-30 day	<0.7-2.5
31-60 day male	0.8-5.0
31-60 day female	0.5-2.3
Children	
3-14 years	0.07-1.7
Females:	
Follicular phase	0.2-1.3
Luteal phase	1.0-4.5
0.2-0.9	
Postmenopausal	

16. REFERENCES

1. Yeo, James K.H. et al. An Automated Solid-Phase 17α -hydroxyprogesterone ELISA method using a microtiter plate. *Ann Clin Biochem* 1988;25:293-297.
2. Hofman, L.F., et al. Direct Solid-Phase Radioimmunoassay for Screening 17α -Hydroxyprogesterone in Whole-Blood Samples from Newborns. *Clin Chem.* 1985;31/7:1127.1130.
3. Sippel W.G., et al. Plasma Levels of 17α -hydroxyprogesterone, Cortisol, and Cortisone During Infancy and Childhood. *Pediatr Res* 1980;14:39-46.
4. Thorneycroft, I.H., et al. The Relation of Serum 17α -hydroxyprogesterone and Estradiol 17B During the Human Menstrual Cycle. *Am J Obstet Gynecol* 1971;111:947-951.

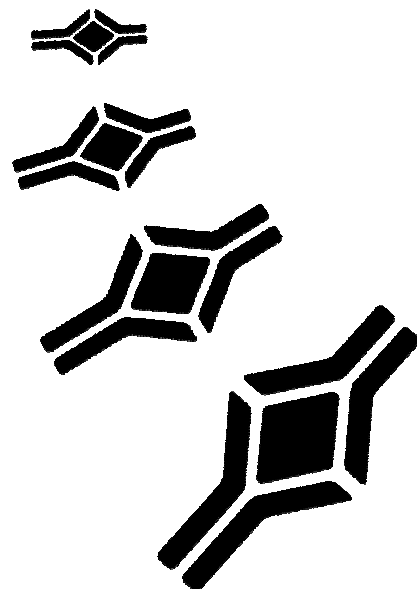

HEADQUARTERS:
 BioVendor Laboratorní
 medicína, a.s.

CTPark Modrice
 Evropska 873

664 42 Modrice
 CZECH REPUBLIC

Phone: +420-549-124-185
 Fax: +420-549-211-460

E-mail: info@biovendor.com
 Web: www.biovendor.com

EUROPEAN UNION:
 BioVendor GmbH

Im Neuenheimer Feld 583

D-69120 Heidelberg
 GERMANY

Phone: +49-6221-433-9100
 Fax: +49-6221-433-9111

E-mail: infoEU@biovendor.com

USA, CANADA AND MEXICO:
 BioVendor LLC

1463 Sand Hill Road
 Suite 227

Candler, NC 28715
 USA

Phone: +1-828-670-7807
 +1-800-404-7807
 Fax: +1-828-670-7809

E-mail: infoUSA@biovendor.com

CHINA - Hong Kong Office:
 BioVendor Laboratories Ltd

Room 4008
 Hong Kong Plaza, No.188

Connaught Road West
 Hong Kong, CHINA

Phone: +852-2803-0523
 Fax: +852-2803-0525

E-mail: infoHK@biovendor.com

CHINA – Mainland Office:
 BioVendor Laboratories Ltd

Room 2405 YiYa Tower
 TianYu Garden, No.150

Lihe Zhong Road
 Guang Zhou, CHINA

Phone: +86-20-87063029
 Fax: +86-20-87063016

E-mail: infoCN@biovendor.com