



AssayMax Human Urokinase (uPA) ELISA Kit

Catalog # EU1001-1

Introduction

Urokinase-type plasminogen activator (uPA) is a highly restricted serine protease that converts the zymogen plasminogen to active plasmin, a broad-spectrum serine proteinase capable of degrading most of the major protein components of the extracellular matrix. Binding of uPA to its receptor and subsequent uPA mediated pericellular proteolysis are involved in many process including cell migration and tissue remodeling in angiogenesis, atherogenesis, tumor cell metastasis, and ovulation (1, 2). High level of uPA is a poor prognostic marker for aggressive breast cancer, aggressive prostate cancer, bladder cancer and gastric cancer (3, 4, 5).

Principal of the Assay

The AssayMax Human Urokinase (uPA) ELISA kit is designed for detection of human uPA in plasma, serum, tissue, cell culture supernatants and urine. This assay employs a quantitative sandwich enzyme immunoassay technique which measures urokinase in 4 hours. A polyclonal antibody specific for uPA has been pre-coated onto a microplate. Urokinase in standards and samples is sandwiched by the immobilized antibody and a biotinylated polyclonal antibody specific for uPA, which is recognized by a streptavidin-peroxidase conjugate. All unbound material is then washed away and a peroxidase enzyme substrate is added. The color development is stopped and the intensity of the color is measured.

Caution and Warning

- This kit is for research use only.
- The kit should not be used beyond the expiration date.
- The Stop Solution is an acid solution.

Reagents

- **uPA Microplate:** A 96-well polystyrene microplate (12 strips of 8 wells) coated with a polyclonal antibody against uPA.
- **Sealing Tapes:** Each kit contains 3 pre-cut, pressure-sensitive sealing tapes that can be cut to fit the format of the individual assay.
- **uPA Standard:** Human uPA in a buffered protein base (24 ng, lyophilized).
- **Biotinylated uPA Antibody (100x):** A 100-fold biotinylated polyclonal antibody against human uPA (80 µl).

- **Streptavidin-Peroxidase Conjugate (SP Conjugate):** A 100-fold concentrate (90 µl)
- **MIX Diluent Concentrate (10x):** A 10-fold buffered protein base (30 ml).
- **Wash Buffer Concentrate (10x):** A 10-fold concentrated buffered surfactant (2 x 30 ml).
- **Chromogen Substrate:** A ready-to-use stabilized peroxidase chromogen substrate tetramethylbenzidine (8 ml).
- **Stop Solution:** A 0.5 N hydroxychloric acid (12 ml) to stop the chromogen substrate reaction.

Storage Condition

- Store unopened kit at 2-8⁰C up to expiration date.
- Opened reagents may be stored for up to 1 month at 2-8⁰C. Store reconstituted standard at -20⁰C or below.
- Opened unused strip wells may return to the foil pouch with the desiccant pack, reseal along zip-seal. May be stored for up to 1 month in a vacuum desiccator.

Other Supplies Required

- Microplate reader capable of measuring absorbance at 450 nm
- Pipettes (1-20 µl, 20-200 µl, and multiple channel)
- Deionized or distilled reagent grade water

Sample Collection, Preparation and Storage

- **Plasma:** Collect plasma using one-tenth volume of 0.1 M sodium citrate as an anticoagulant. Centrifuge samples at 2000 x g for 10 minutes and use supernatants. Dilute samples 1:2 with MIX Diluent and assay. The undiluted samples can be stored at -20⁰C or below for up to 3 months. Avoid repeated freeze-thaw cycles. (EDTA or Heparin can also be used as anticoagulant.)
- **Serum:** Samples should be collected into a serum separator tube. After clot formation, centrifuge samples at 2000 x g for 10 minutes. Remove serum and assay. Dilute samples 1:2 into MIX Diluent. Store serum at -20⁰C or below. Avoid repeated freeze-thaw cycles
- **Cell Culture Supernatants:** Centrifuge cell culture media at 2000 x g for 10 minutes to remove debris. Collect supernatants and assay. The samples can be stored at -20⁰C or below. Avoid repeated freeze-thaw cycles.
- **Tissue:** Extract tissue samples with 0.1 M phosphate-buffered saline (pH7.4) containing 1% Triton X-100 and centrifuge at 14000 x g for 30 min. Collect the supernatant and measure the protein concentration. Dilute the tissue extract 1:2 into MIX Diluent and assay. The undiluted samples can be stored at -20⁰C or below for up to 3 months. Avoid repeated freeze-thaw cycles.
- **Urine:** Collect urine using sample pot. Centrifuge samples at 800 x g for 10 minutes and use supernatants. Dilute samples 1:20 into MIX Diluent and assay. The undiluted samples can be stored at -20⁰C or below for up to 3 months. Avoid repeated freeze-thaw cycles.

Reagent Preparation

- Freshly dilute all reagents and bring all reagents to room temperature before use. If crystals have formed in the concentrate, mix gently until the crystals have completely dissolved.
- **Standard Curve:** Reconstitute the 24 ng of human uPA Standard with 3 ml of MIX Diluent to generate an 8 ng/ml of stock solution. Allow the standard to sit for 10 minutes with gentle agitation prior to making dilutions. Prepare triplicate standard points by serially diluting the uPA standard solution (8 ng/ml) twofold with equal volume of MIX Diluent to produce 4, 2, 1, 0.5, 0.25, and 0.125 ng/ml solutions. MIX Diluent serves as the zero standard (0 ng/ml).

Standard Point	Dilution	[uPA] (ng/ml)
P1	1 part Standard (8 ng/ml)	8.000
P2	1 part P1 + 1 part MIX Diluent	4.000
P3	1 part P2 + 1 part MIX Diluent	2.000
P4	1 part P3 + 1 part MIX Diluent	1.000
P5	1 part P4 + 1 part MIX Diluent	0.500
P6	1 part P5 + 1 part MIX Diluent	0.250
P7	1 part P6 + 1 part MIX Diluent	0.125
P8	MIX Diluent	0.000

- **MIX Diluent Concentrate (10x):** Dilute MIX Diluent Conc. 1:10 with reagent grade water.
- **Biotinylated uPA Antibody (100x):** Spin down the antibody briefly and dilute the desired amount of the antibody 1:100 with MIX Diluent.
- **Wash Buffer Concentrate (10x):** Dilute Wash Buffer Conc. 1:10 with reagent grade water.
- **SP Conjugate (100x):** Spin down the SP Conjugate briefly and dilute the desired amount of the conjugate 1:100 with MIX Diluent.

Assay Procedure

- Prepare all reagents, working standards and samples as instructed. Bring all reagents to room temperature before use. The assay is performed at room temperature (20-30°C).
- Remove excess microplate strips from the plate frame and return them immediately to the foil pouch with desiccant inside. Reseal the pouch securely to minimize exposure to water vapor and store in a vacuum desiccator.
- Add 50 µl of Standard or sample per well. Cover wells and incubate for two hours. Start the timer after the last sample addition.
- Wash five times with 200 µl of Wash Buffer. Invert the plate and decant the contents, and hit it 4-5 times on absorbent paper towel to completely remove liquid at each step.
- Add 50 µl of Biotinylated uPA Antibody to each well and incubate for 60 minutes.
- Wash five times with 200 µl of Wash Buffer as above.
- Add 50 µl of Streptavidin-Peroxidase Conjugate per well and incubate for 30 minutes. Turn on the microplate reader and set up the program in advance.
- Wash five times with 200 µl of Wash Buffer as above.
- Add 50 µl of Chromogen Substrate per well and incubate for about 10 minutes or till the optimal blue color density develops. Gently tap the plate to ensure thorough mixing and break the bubbles in the well with pipette tip.
- Add 50 µl of Stop Solution to each well. The color will change from blue to yellow.


- Read the absorbance on a microplate reader at a wavelength of 450 nm **immediately**. Please note that some unstable black particles may be generated at high concentration points after stopping the reaction for about 10 minutes, which will reduce the readings.

Data Analysis

- Calculate the mean value of the triplicate readings for each standard and sample.
- To generate a standard curve, plot the graph using the standard concentrations on the x-axis and the corresponding mean 450 nm absorbance on the y-axis. The best-fit line can be determined by regression analysis using 4-parameter or log-log curve fit.
- Determine the unknown sample concentration from the Standard Curve and multiply the plasma or tissue mean value by the dilution factor.

Standard Curve

- The curve is provided for illustration only. A standard curve should be generated each time the assay is performed.

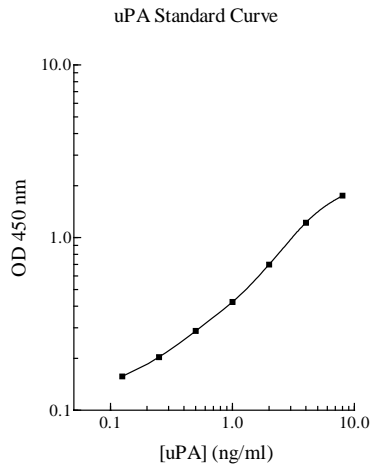

Precision, Sensitivity and Specificity

- The minimum detectable level of uPA is typically < 100 pg/ml.
- Intra-assay and inter-assay coefficients of variation were 5.2 % and 7.5% respectively.
- This assay recognizes single chain, two-chain, and both receptor and PAI-bound human uPA.

Cross-Reactivities

Species	% Cross Reactivity
Monkey	< 20 (suggest dilution 1:2 for plasma)
Mouse	< 10
Rat	None
Swine	< 20 (suggest dilution 1:2 for plasma)
Beagle	< 10
Bovine	None

Recovery

Standard Added Value	0.2 – 2 ng/ml
Recovery %	87-117 %
Average Recovery %	102 %

Linearity

Sample Dilution	Average Percentage of Expected Value	
	Plasma	Serum
1:2	100%	100%
1:4	98%	103%
1:8	101%	100%

Sample Dilution	Average Percentage of Expected Value
	Urine
1:10	100%
1:20	97%
1:40	101%

References

1. Okada, S. *et al.* (1996) *Arterioscl. Thromb. Vasc. Biol.* 16: 1269
2. Besser, D. *et al.* (1996) *Fibrinolysis* 10: 215
3. Duffy, M.J. *et al.* (1990) *Cancer Res.* 50:6827
4. Hasui, Y. *et al.* (1992) *Int. J. Cancer* 50: 871
5. Nishino, N. *et al.* (1988) *Thromb. Res.* 50:527

Revision 7.1