

**STAPHYLOCOCCUS AUREUS
(mecA-) DNA CONTROL**

FOR INFORMATION USE ONLY

Not to be used for performing the assay. Refer to the insert accompanying the kit.

MBC073: Purified DNA of *Staphylococcus aureus* (mecA-).

LOT SPECIFICATIONS:

Microorganism: *Staphylococcus aureus* (mecA-).

Preparation: Grown in blood agar culture medium.

TECHNICAL DATA:

Presentation: Lyophilized.

Concentration: 50 ng/μl when reconstituted.

Instructions for reconstitution:

- 1.-Centrifuge 1 minute at 1000 g.
- 2.-Add 50 μl of sterile bidistilled water and mix until completely resuspended. The concentration will be 50 ng/μl.
- 3.-Shake with vortex for 30 seconds to dissolve and homogenize completely.

Extract preparation: SDS and proteinase K treatment, phenol/chloroform extraction and ethanol precipitation.

Shipping and storage:

Special transport conditions not required. Reconstitute upon receipt. Store between -5°C and -40°C after reconstitution. Avoid multiple freeze-thaw cycles as product degradation may result.


Recommendations and precautions:

This product is not for diagnostic use.

This product is for research use only.

This lot of DNA was prepared from *Staphylococcus aureus* (mecA-). The material should be handled and disposed of as potentially infectious. Observe the local regulations for clinical waste disposal.

PCR analysis of DNA control: PCR analysis was performed with a specific oligo pair on 50 ng of purified *Staphylococcus aureus* (mecA-) DNA. The reaction produced a 94 bp fragment. It was visualized on a 2% agarose gel using ethidium bromide staining. The gel photograph is shown below:

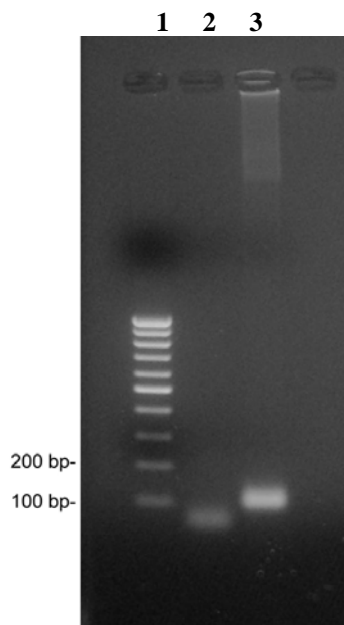

Line	Sample
1	Molecular size marker (100 bp)
2	Negative control
3	94 bp amplified fragment

For inquiries please contact:
customerservice@vircell.com

REVISED: 07/2008

FOR RESEARCH USE ONLY

Manufacturer: VIRCELL, S.L. Pza. Dominguez Ortiz 1. Polígono Industrial Dos de Octubre.18320 Santa Fe *GRANADA* SPAIN* Tel.+34.958.441264* Fax+34.958.510712
<http://www.vircell.com>